



THUASNE

ROM-R



Unique size

U04602

Target patient

The brace is designed for patients with severe ligament injuries, or for post-operative treatments.

Indications

- Post-operative and/or rehabilitative knee immobilization (including after ligament injuries, patellar tendon tears, meniscus tears, cartilage injuries, as well as tibia plateau and condylar fractures)
- Post-traumatic knee immobilization
- Conservative treatment of the knee after ligament tears / injuries
- Knee instability

ROM-R is a regulated medical device which bears CE marking. Read all instructions before use and seek advice from a qualified healthcare professional.

Contraindications: please refer to the instructions for use.



Thuasne - 120, rue Marius Auphan
92300 Levallois Perret, France



© Studio Caterin - Jennifer Williams - 2001009 - 2020-01

ROM-R

POST-OPERATIVE

KNEE SPLINT

WITH A TOOL-LESS

ADJUSTMENT

SYSTEM




ROM-R

Post-operative knee splint with a tool-less adjustment system


EASE OF USE

- Tool-less length adjustment system for a better fitting (proximal and distal)

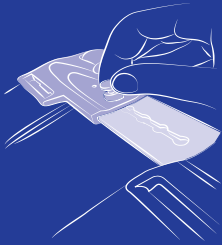


- 4 quick buckles to put the brace on and off

Lock

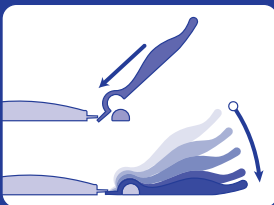

Insert vertically and flip down

Unlock

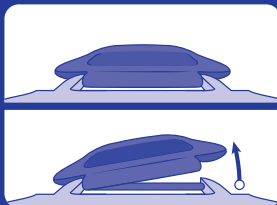

Twist buckle up from one corner


COMFORT

- Foam on surfaces in direct contact with the skin

MODE OF ACTION

- Immobilization provided by 2 lateral uprights

- **Drop Lock system to immobilize the leg** by 15° steps from 0° to 90°

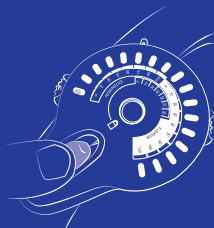

- Flexion and extension can be settled without any tool by 15° steps in flexion (from 0° to 120°) and extension (from 0° to 105°)

